

The OZOMED smartline is characterized by its highly compact design and excellent performance. With an adjustable ozone concentration between 1 and  $\leq 80 \mu\text{g/ml}$ , all common ozone treatments can be performed. The ozonization of water is possible with the optionally available water ozonization unit.

## Product highlights

Easy handling via colored touchscreen, supporting multiple language options.

As a tabletop unit, ideal for applications in small and medium-sized medical practices.

Powerful catalyst for safe ozone destruction.

Low maintenance with above-average product service life.

Developed and manufactured in Germany.



Technical specification	OZOMED smartline
Device setup	Tabletop unit
Mode of operation	Normobaric infusion technique
Ozone concentration	1 - 80 $\mu\text{g/ml}$
Flow rate	0.9 l/min
Operating pressure	1 bar
Vacuum pump	Up to -0.7 bar
Water treatment with ozone	Optionally available
Dimensions [W x H x D]	42 x 19 x 39 cm
Weight	9.8 kg
Electrical data	100 - 240 V AC; 50 - 60 Hz; 1.2 A; Protection class I
Classification	IIb medical device according to EU Directive 93/42/EEC
Warranty	2 years in Europe

Form of therapy	OZOMED smartline
Major autohemotherapy	●
Minor autohemotherapy	●
Rectal insufflation	●
Bag treatment	●
Low pressure ozone gas application	
Ozone injections	●
Ozonized water	optionally




## Convenience kits for ozone therapy

### OZONOSAN Set 4 Plus Convenience kit



- 1 OZONOSAN vacuum bottle 250 ml, microbubble system, 12 ml sodium citrate as anticoagulant, latex-free
- 1 Blood transfusion set with steel tip
- 1 Perfusion set optionally with 0.8 or 1.1 mm needle
- 1 Disposable syringe 50 ml, siliconized
- 1 Transfer filter set with cannula and filter (0.2 µm) for ozone introduction

### OZONOSAN Set Rectal insufflation with dosing bellow

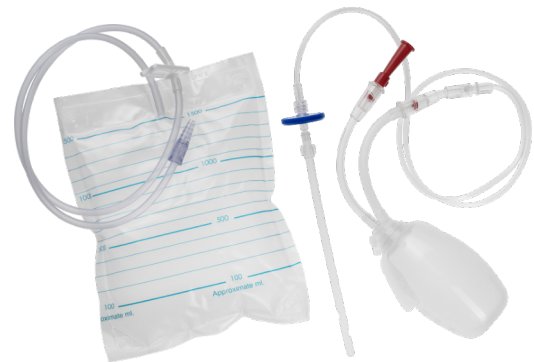

- 1 Dosing bellow
- 5 Ozone reservoirs with valve
- 50 Special disposable catheters 18 Ch.
- 50 Infusion lines 75 cm
- 50 Bacteria filter

## Hänsler Medical GmbH

Nordring 8, 76473 Iffezheim, Germany

+49 (0) 7229 3046 0

@ info@haensler-medical.com

www.haensler-medical.com

