

OZONOSAN GAMMA

DIE NEUESTE ENTWICKLUNG VON HÄNSLER MEDICAL FÜR ALLE THERAPIEFORMEN


THE LATEST DEVELOPMENT FROM HÄNSLER MEDICAL FOR ALL FORMS OF OZONE THERAPY

Das neueste Ozontherapiesystem von Hänsler Medical kombiniert bewährte OZONOSAN-Qualität mit elegantem Design und Funktionalität. Einstellbare Ozonkonzentrationen von 1 bis $\leq 80 \mu\text{g/ml}$ ermöglichen alle gängigen Anwendungsformen – von Großer Eigenblutbehandlung, Injektionsverfahren, Beutelbegasung bis hin zu Wasserozonisierung mit der optional erhältlichen Wasserozonisierungseinheit.

The OZONOSAN gamma is the latest product innovation from Hänsler Medical. It combines proven OZONOSAN quality with remarkable design and functionality. With an adjustable ozone concentration between 1 and $\leq 80 \mu\text{g/ml}$, the OZONOSAN gamma enables all common types of ozone therapy at your fingertips via large easy to navigate touchscreen.

HIGHLIGHTS

- Einfache Bedienung über farbigen Touchscreen, mehrere Sprachoptionen werden unterstützt. / *Easy handling via colored touchscreen, supporting multiple language options.*
- Als Tischgerät ideal für Anwendungen in kleinen und mittelgroßen Arztpraxen. / *As a tabletop unit, ideal for applications in small and medium-sized medical practices.*
- Wartungsarm bei gleichzeitig hoher Produktlebensdauer. / *Low maintenance with a long product life.*
- Premium Design für die gehobene Praxis. / *Premium design for the professional clinic.*
- In Deutschland entwickelt und produziert. / *Developed and manufactured in Germany.*

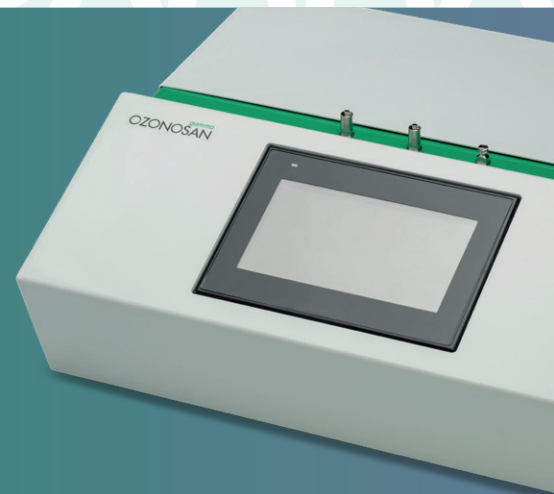

Technische Spezifikation / Technical Specifications

Aufbau / Device setup

Eigenbluttherapie / *Major Autohemotherapy type*

Konzentrationsbereich / Ozone concentration range

Volumenstrom / Flow rate

Vakuumpumpe / Vacuum pump

Wasseraufbereitung mit Ozon / Water ozonization

Abmessung [H x B x T] / Dimensions [H x W x D]

Gewicht / Weight

Elektrische Daten / Electrical data

EU Medizinprodukte Zertifizierung / EU Medical Device Certification

Gewährleistung / Warranty

OZONOSAN gamma

Tischgerät, kompakt / *Tabletop unit*

normobare Infusionstechnik im offenen oder geschlossenen System / *normobaric infusion technique in an open or closed system*

1 - 80 $\mu\text{g/ml}$ / *1 - 80 $\mu\text{g/ml}$*

0,9 l/min / *0.9 l/min*

bis zu -0,7 bar / *Up to -0.7 bar*

OZONOSAN Wasserozonisierungseinheit optional erhältlich / *OZONOSAN water ozonization unit optionally available*

18 x 47 x 46 cm / *18 x 47 x 46 cm*

10,7 kg / *10,7 kg*

U: 100 - 240 V AC, P: 200 - 200 VA, f: 50 - 60 Hz; Schutzklasse I / *U: 100 - 240 V AC, P: 200 - 200 VA, f: 50 - 60 Hz; Protection class I*


MDR (CE0483) / *MDR (CE0483)*

2 Jahre / *2 years*

OZONOSAN^{gamma}

Anwendungsformen / Form of therapy

- ✓ Große Eigenblutbehandlung / *Major autohemotherapy*
- ✓ Rektale Insufflation / *Rectal insufflation*
- ✓ Beutelbegasung / *Bag gassing*
- ✓ Ozon-Injektionen / *Ozone injections*
- ! Ozonisiertes Wasser (optional) / *Ozonized water (optionally)*

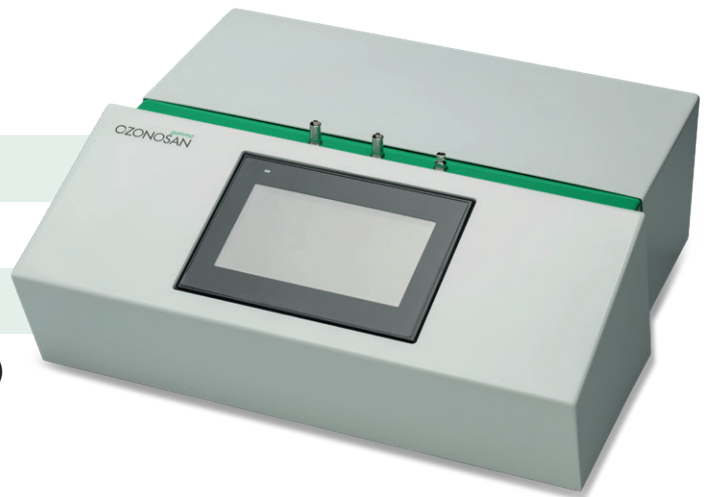

Convenience Kits für die Ozontherapie / Convenience kits for ozone therapy



OZONOSAN Set 4 Plus Convenience Kit

OZONOSAN Set 4 Plus Convenience kit

- 1 OZONOSAN Vakuumflasche 250 ml, Mikroperl-System, 12 ml Natriumcitrat als Gerinnungshemmer, latexfrei /
- 1 OZONOSAN vacuum bottle 250 ml, microbubble system, 12 ml sodium citrate as anticoagulant, latex-free
- 1 Transfusionsbesteck mit Stahlspitze /
- 1 Blood transfusion set with steel tip
- 1 Perfusionsbesteck wahlweise 0,8 oder 1,1 mm Nadel /
- 1 Perfusion set optionally with 0.8 or 1.1 mm needle
- 1 Einmalspritze 50 ml silikonisiert /
- 1 Disposable syringe 50 ml, siliconized
- 1 Transferfilterset /
- 1 Transfer filter set with cannula and filter (0.2 µm) for ozone introduction



Sangiset Plus Convenience Kit

Sangiset Plus Convenience Kit

- 1 Sangiset, 300ml /
- 1 Sangiset, 300ml
- 1 Ampulle Natrium Citrat steril, 10ml /
- 1 Ampoule sodium citrate, sterile, 10ml
- 1 Perfusionsbesteck mit Präzisionsschliff, 1,1mm /
- 1 Perfusions set (butterfly) with precision cutting, 1.1mm